(1) UBIGER

(1) UBIGENE **STR Authentication Report For Cell Line**

Tel: +86 153 6067 3248(Int'l) US Toll free: 855 777 3210 EU Toll free: 800 3272 9252

Korea Toll free: 001 800 3272 9252

Email: info@ubigene.com Web: www.ubigene.us

Online Inquiry

1. Cell Name: Detroit 562

2. Test Method: DNA was extracted using the genome extraction kit (Axygen), amplified using a 21- STR amplification protocol, the STR loci and gender gene Amelogenin were detected on the ABI 3730XL Genetic Analyzer.

| , | Genotype | analysis resu | Its of STR and | l Amelogenin l | oci in cells | | |
|----------|---|---------------|----------------|--|--------------|--------|--|
| STR Loci | STR profile for Sample cell Cell name: Detroit 562 | | | STR profile for Cell Bank cell Cell name: Detroit-562 | | | |
| | | | | | | | |
| | D5S818 | 11 | 12 | UBIC | 11 | 12 | |
| D13S317 | 12 | 12 | | 12 | 12 | -< | |
| D7S820 | 8 | 10 | | 8 | 10 | -11BIG | |
| D16S539 | 11 | 11 | | 11 | 11 | | |
| VWA | 16 | 16 | | 16 | 16 | | |
| TH01 | 8 | 9 | | 8 | 9 | | |
| AMEL | Х | Х | | Х | Х | | |
| ТРОХ | 8 | 10 | .0 | 8 | 10 | | |
| CSF1PO | 11 | 13 | UBIL | 11 | 13 | | |
| D12S391 | 17 | 18 | | | | ~~ | |
| FGA | 21 | 21 | | | | · IBIO | |
| D2S1338 | 25 | 25 | | | | | |
| D21S11 | 28 | 30 | | | | | |
| D18S51 | 15 | 15 | | | | | |
| D8S1179 | 13 | 14 | | | | | |

Tel: +86 153 6067 3248(Int'l) US Toll free: 855 777 3210 EU Toll free: 800 3272 9252

Korea Toll free: 001 800 3272 9252

Email: info@ubigene.com Web: www.ubigene.us Online Inquiry

| D3S1358 | 15 | 16 | | | |
|---------|------|------|------|-----|--------|
| D6S1043 | 12 | 18 | | | |
| PENTAE | 13 | 13 | . C | ENE | |
| D19S433 | 14 | 14 | UBIC | | |
| PENTAD | 13 | 13 | | | 757 |
| D1S1656 | 15.3 | 17.3 | | | * IBIO |


4. Conclusion: The STR results of <u>Detroit 562</u> cells are consistent with the genotypes of <u>Detroit 562</u> cell lines in <u>DSMZ</u> database, the cell ID corresponded to <u>CCL-138</u>, and the STR results <u>completely matched</u>. In the test, <u>No multi allele was found</u> in this cell line.

Tel: +86 153 6067 3248(Int'l) US Toll free: 855 777 3210 EU Toll free: 800 3272 9252

Korea Toll free: 001 800 3272 9252

Email: info@ubigene.com Web: www.ubigene.us Online Inquiry

5. Attached Image

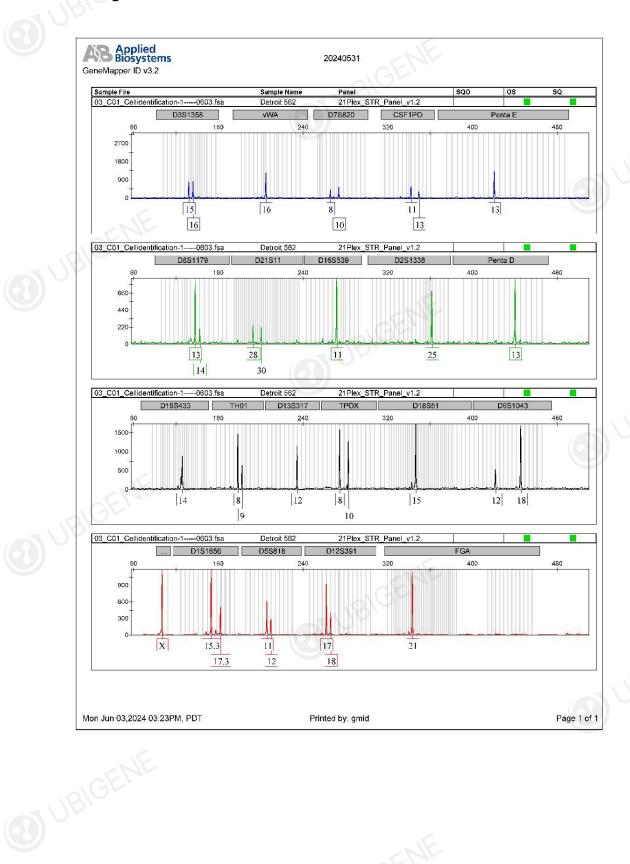

Tel: +86 153 6067 3248(Int'l) US Toll free: 855 777 3210 EU Toll free: 800 3272 9252

Korea Toll free: 001 800 3272 9252

Email: info@ubigene.com
Web: www.ubigene.us

Online Inquiry